



User Guide

CentrePoint Insight Watch

Activity Monitor: CentrePoint Insight Watch | E.200.6002 | Revision: 4 | Released: 02/18/2020



CentrePoint Insight Watch



Activity Monitor: CentrePoint Insight Watch | E.200.6002 | Revision: 4 | Released: 02/18/2020

Contents

Overview.....	2
What You'll Need.....	3
Getting Started.....	4
Assigning the Insight Watch to a Subject.....	6
Information for Subjects.....	11
Uploading Data to CentrePoint.....	13
Unassigning the Insight Watch.....	14
Lost or Broken Activity Monitors.....	15
Frequently Asked Questions.....	17
Regulatory Information.....	18

Overview



CentrePoint Insight Watch

The CentrePoint (CP) Insight Watch from ActiGraph captures and records continuous, high resolution acceleration data, which are converted into variety of objective physical activity and sleep measures within the cloud-based CentrePoint software system. This manual provides instructions on how to setup, deploy, and upload data from the CentrePoint Insight Watch to the CentrePoint software system.

Specifications

Dimensions	5.01 x 3.46 x 1.05 cm	Dynamic range	+/- 8G
Weight	14 grams	Communication	USB, Bluetooth® LE 5
Sample rate	32 Hertz (default) - 256 Hertz	Water resistance	1 meter, 30 minutes
Battery life	30 days*	Wear location	Wrist
Data storage	30 days/512 MB	Warranty	1 year

*At 32 Hz (default) sample rate.

What You'll Need



PC + ActiSync Software



CentrePoint Account



CentrePoint Insight Watch



CentrePoint Data Hub + Power Cable (optional)



Charging Dock



USB Cable



Wall Charger (optional)

1

Allow firewall access to the following URLs

Contact your local IT department prior to assigning activity monitors to ensure you have firewall access to the URLs listed below.

ActiSync service:

<https://studyadmin-actisync-service.actigraphcorp.com>

Main Web Portal:

<https://centrepont.actigraphcorp.com>

Blob Azure Storage:

<https://acticlouduploadsprod.blob.core.windows.net>

2

CentrePoint account setup

Using the invitation link provided by ActiGraph or your study sponsor, create a CentrePoint username and password.



3

Install ActiSync software

- 1 Download the ActiSync software from:
<https://actigraphcorp.com/actisync>
- 2 Double click the file to open and follow prompts to complete installation.



System Requirements

OS Windows 7 or later | **CPU** Dual core processor | **HDD** 20MB, USB 2.0

4

Charge Insight Watch

The Insight Watch should be fully charged before assigning it to a subject. It takes 2-3 hours to fully charge the device.

- 1 Connect the charging dock to the PC or a wall outlet using the appropriate cable.
- 2 Connect the Insight Watch to the charging dock by pulling the plastic tab away from the dock and inserting the device with the screen facing up and the button on the opposite side of the cable port.
- 3 The Insight Watch is fully charged once a solid battery icon appears on the screen.



Battery Life Specifications: The battery life of the Insight Watch is approximately 30 days.*

Note: An Internet connection is not required to charge the activity monitors.

*At 32 Hz (default) sample rate.

Assigning the Insight Watch to a Subject


1

Assign the Insight Watch to a new subject

- 1 Verify that ActiSync is running on the PC by checking for a small icon in the taskbar.



- 2 Connect the charging dock to the computer using the provided USB cable and insert the Insight Watch into the charging dock.
- 3 Once connected, the ActiSync icon will turn yellow and a dialog box will open. Select “Assign the subject now.”

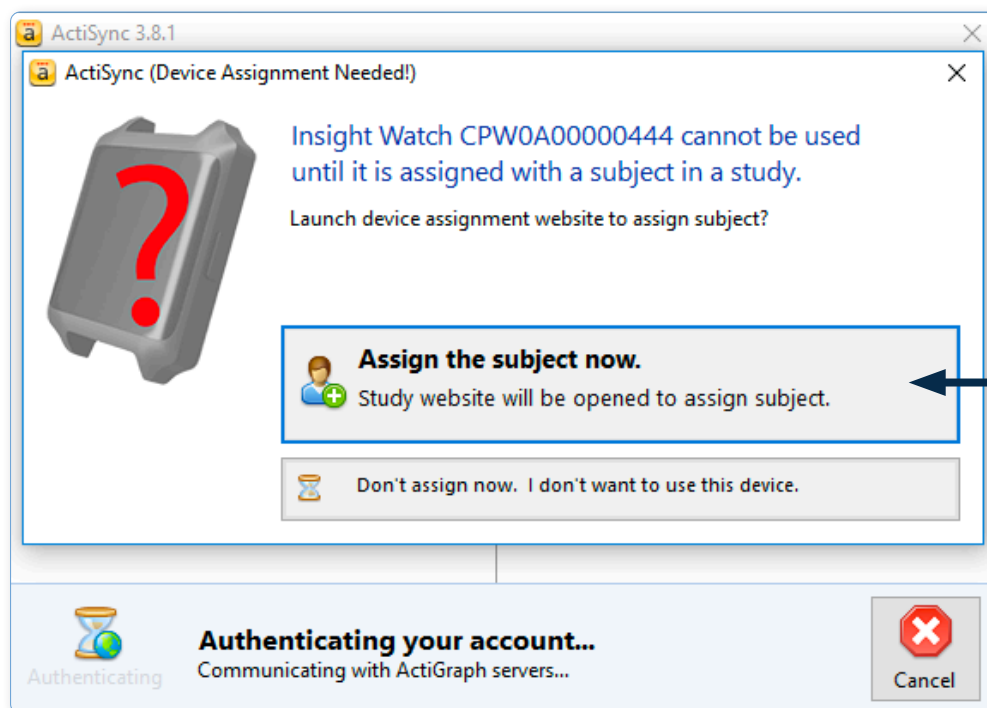

- 4 Enter your CentrePoint login information to open the CentrePoint web page in your browser.

Assigning the Insight Watch to a Subject (Continued)

- 5 Complete the **Device Assignment** form and select “Assign Device to a New Subject.” Wait for progress bar to indicate that assignment is complete.

The screenshot shows the CENTREPOINT interface for device assignment. On the left is a sidebar with navigation options: Choose A Study, Dashboard, Manage Site Users, Milestones, Subjects, Devices (highlighted), Logistics, Exports, Help, Internal Systems, and Administration. The main content area is titled 'Device Assignment' and shows a battery status for device CPW1A37180213 at 99% (4.13V). Below this is a section '+ Assign To a New Subject' with a form. The form includes: 'Selected Study' (Demo Organization/Demo Study), 'Choose Site from Study' (Demo Study / Site 1), 'Subject Identifier' (101), 'DOB' (02/09/1980), 'Weight' (155.0 lbs), and 'Wear Position' (Left Non-Dominant Wrist). A blue arrow points to the 'Assign Device to a New Subject' button.

Device Assignment

99% The battery of the connected activity monitor, CPW1A37180213 is 99% (4.13V)
Monitor last seen a few seconds ago.

+ Assign To a New Subject

Selected Study Demo Organization/Demo Study

Choose Site from Study Demo Study / Site 1

Subject Identifier 101

DOB 02/09/1980

Weight 155.0 lbs

Wear Position Left Non-Dominant Wrist

Assign Device to a New Subject

The screenshot shows the 'Activity Monitor Assignment' completion screen for device CPW0A00000444. It displays a 99% battery status (4.16V). The 'Assignment Complete!' section includes 'Next Steps' (activity monitor ready for Subject 1234AN) and a note about assigning another subject. A three-step process diagram shows: 1. Activity monitor connected to a computer, 2. Waiting (hourglass icon), and 3. Activity monitor on a hand. At the bottom, a green bar indicates 'Activity Monitor Configured' and provides details: Serial Number CPW0A00000444, Subject 001002424, and Data Hub CDM1A11180099.

Activity Monitor Assignment Device: CPW0A00000444

99% The battery of the connected activity monitor, CPW0A00000444 is 99% (4.16V)
Monitor last seen a few seconds ago.

Assignment Complete!

Next Steps

- The activity monitor is ready to be placed on [Subject 1234AN](#) for data collection.

Need to assign another subject?

Simply close this window, or tab, and connect your next activity monitor to ActiSync. Follow the prompt to assign the activity monitor to a subject.

1 2 3

Activity Monitor Configured

You may unplug the activity monitor. The activity monitor can now be placed on the subject.


Serial Number CPW0A00000444

Subject 001002424

Data Hub CDM1A11180099

Assigning the Insight Watch to a Subject (Continued)

Note: To assign an activity monitor to an existing subject (such as in the case of a lost or broken device), complete steps 1-4 above, and scroll to the bottom of the page. Select the corresponding subject identifier from the list, then select **"Assign Device to an Existing Subject."** Wait for the progress bar to indicate that assignment is complete.

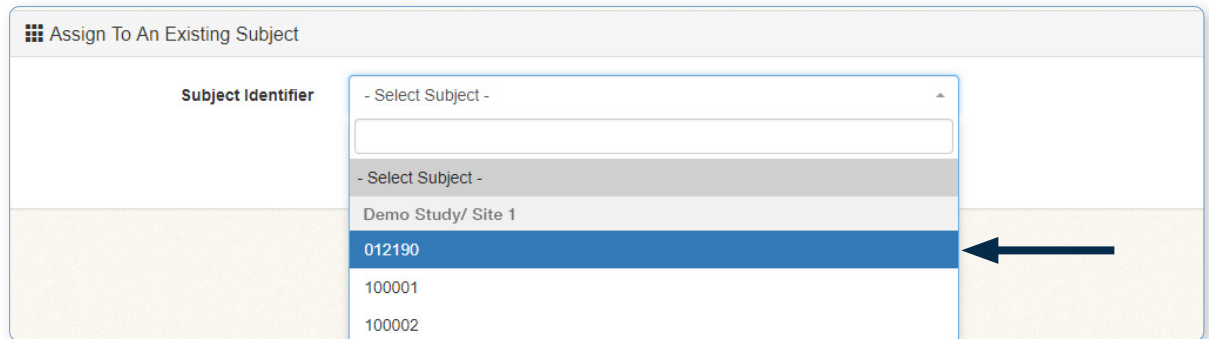

Assign To An Existing Subject

Subject Identifier

- Select Subject -


- Select Subject -

Demo Study/ Site 1

012190

100001

100002

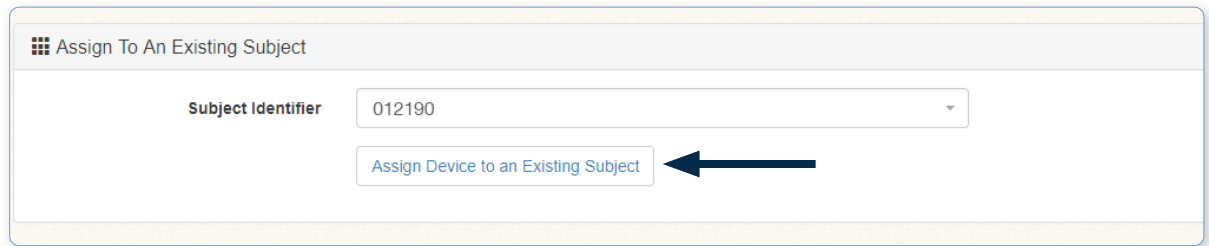

Assign To An Existing Subject

Subject Identifier

012190

Assign Device to an Existing Subject

- 6 Remove the Insight Watch from the charging dock. Wait for the activity monitor to go from **halt mode** into **active mode**, as indicated by the running man icon in the upper left corner of the screen. This takes approximately 90 seconds. The Insight Watch is now collecting data and can be given to the subject.



Halt

Stop icon will be displayed if the Insight Watch is not collecting data.



Active Mode

The Insight Watch is collecting data.

2

Pair the Insight Watch with a CentrePoint (CP) Data Hub (If you are not using a CP Data Hub, skip ahead to **Information for Subjects**)

The CP Data Hub passively acquires collected raw accelerometry data from the Insight Watch and transmits to the cloud-based CentrePoint system over a cellular network.

- 1 Connect the adapter end of the power cable into the back of the CP Data Hub and plug the other end into a wall outlet.



- 2 Connect the charging dock to the USB port on the CP Data Hub using the USB cable.
- 3 Insert the Insight Watch into the charging dock. Once connected, a progress wheel will appear on the screen indicating that an upload has begun. When the upload is complete, a checkmark will appear on the screen.



Upload in Progress



Upload Complete



Warning Sign

If you see a warning sign on the screen of the Insight Watch, remove it from the charging dock and insert again. Make sure you see a checkmark on the screen before proceeding to the next step.

Assigning the Insight Watch to a Subject (Continued)

- 4 Verify connection by logging into your CentrePoint account and selecting the “Subjects” tab from the left menu. Select the desired subject from the list.

The screenshot shows the 'Study Subjects' page in the CentrePoint interface. The left sidebar has the 'Subjects' tab selected. The main content area shows a table with columns: Subject, Site, and Days Of Data. The first subject, 0101700, is highlighted with an arrow. The table lists three subjects:

Subject	Site	Days Of Data
0101700 Female 6/29/1971		20 days
0101800 Female 9/9/1959	Site 01	35 days
01101 Female 2/10/1980	Site 01	12 days

- 5 You will be directed to the “Subject Details” page. Select the “Subject Timeline” tab, and verify that an upload has occurred via Data Hub USB. The Insight Watch and the CP Data Hub are now paired.

The screenshot shows the 'Subject Details' page for subject 111001. The left sidebar has the 'Subjects' tab selected. The main content area shows the 'Subject Timeline' tab. The timeline lists events with UTC Timestamps and descriptions. An arrow points to the 'upload' event.

UTC Timestamp	Event
1/13/2017 7:37:50 PM	upload Upload containing 1 file. Epoch summary file contributed 3 minutes of new activity data. - Battery: 90% (3.88V) - Source: Data Hub USB
1/13/2017 7:34:14 PM	assignment Assignment of device CPW0A00000444 completed by initialization
1/13/2017 7:34:07 PM	assignment Device CPW0A00000444 assigned to subject
1/13/2017 7:34:05 PM	added via assignment StudyAdmin - name.lastname@company.org
	weight entry Subject weight recorded - 150 lbs

- 6 Remove the Insight Watch from the charging dock. Give the Insight Watch, charging dock, CP Data Hub, USB cable, and power cable to the subject.

1

Wearing the CentrePoint Insight Watch

The Insight Watch is typically worn on the non-dominant wrist. The appropriate wear location will be outlined in the study protocol.



- 1 The subject should be instructed to wear the Insight Watch strapped securely to the wrist with the display text facing right-side up when viewed like a wrist watch.
- 2 The backlight on the display is activated by pressing the button on the side of the Insight Watch.
- 3 The Insight Watch can be worn for 24 hours a day, dependent on the study parameters. There is no need to remove the Insight Watch during regular bathing or swimming activities.

2

Charging the Insight Watch at home

The Insight Watch has a battery life of approximately 30 days when using the default sample rate of 32 Hz. For data collection periods lasting longer than 30 days, subjects should be instructed to charge the Insight Watch at home.

- 1 Provide the subject with a charging dock and a wall charger.
- 2 Instruct the subject to fully charge the Insight Watch every 3-4 weeks for the duration of the study.

3

Setting up the CP Data Hub at home

(If you are not using a CP Data Hub, skip ahead to **Uploading Data to CentrePoint**)

The CP Data Hub should be placed at home in a location where the subject spends most of his or her time.

- 1 Connect the adapter end of the power cable into the back of the CP Data Hub and plug the other end into a wall outlet.



- 2 The blue LED on the back of the CP Data Hub will begin to flash during startup. Once the CP Data Hub connects to the cellular network, the blue LED will stay on solid.



- 3 Once the CP Data Hub is connected, use the USB cable to connect the charging dock to the CP Data Hub.
- 4 Insert the Insight Watch into the dock for charging or docked data uploads, as determined by the study parameters.

Note: The Insight Watch automatically transmits raw data via Bluetooth® to the CP Data Hub each day. It is not necessary for the subject to dock the Insight Watch to upload data from home.

Uploading Data to CentrePoint

- 1 Verify that ActiSync is running and insert the Insight Watch into the charging dock connected to the PC. Data will automatically upload to the CentrePoint system.
- 2 Once the upload is complete, a browser window will open to display the “Subject Upload Compliance” page. Mark the corresponding milestone and review wear compliance.

The image displays two screenshots of the CentrePoint web application interface, specifically the 'Subject Upload Compliance' page for a subject named 003000009.

Top Screenshot: 'Please Mark a Milestone'

- The page header shows 'Demo Organization/Demo Study' and a 'Change Study' link.
- The subject information is displayed: Subject: 003000009, DOB: 2/9/1980, Gender: Male, Weight: 162 lbs, with a 'Subject Detail Views' dropdown.
- A progress bar at the top indicates the workflow: 'UPLOAD DATA' (completed, green), 'MARK MILESTONE' (current step), and 'VERIFY COMPLIANCE'.
- The main section is titled 'Please Mark a Milestone'.
- Under 'Recommended Milestone', there is a list with 'Visit 1' selected.
- Two large blue arrows point from the 'Visit 1' list item to the 'Mark Milestone' button.
- Below the 'Mark Milestone' button is a link: 'I don't want to mark a Milestone'.
- A 'Mark Different Milestone' link is located at the bottom right.

Bottom Screenshot: 'Marked Milestone' and 'Subject Wear Compliance'


- The 'Marked Milestone' section shows 'Visit 1 - Marked by name.lastname@email.com'.
- The 'Subject Wear Compliance' section features a bar chart titled 'Compliance on Range' showing compliance levels from Mar 02 to Mar 08. The chart shows green bars for Mar 03 and Mar 06, and blue bars for Mar 04, Mar 05, Mar 07, and Mar 08. The x-axis is labeled with dates (10/1/2016 - 10/15/2016) and a note: 'Minimum Daily Wear Compliance = 20 hours'.
- Next to the chart is a circular gauge showing '2 of 7 Days'.
- To the right is a 'Total Compliance' section with a circular gauge showing '193 of 200 Days' (10/1/2016 - 10/7/2016).
- At the bottom, there is a blue button labeled 'I Acknowledge Subject 003000009's Wear Compliance', with a large blue arrow pointing to it from the right.

- 3 Remove the Insight Watch from the charging dock and return it to the subject for continued monitoring.


Unassigning the Insight Watch

Once the data collection period is over, you will need to unassign the Insight Watch from the subject.

- 1 Follow steps 1-2 above to perform a final data upload.
- 2 From the “Subject Details” page, select “Stop Collecting Data” to remove the subject assignment from the Insight Watch.

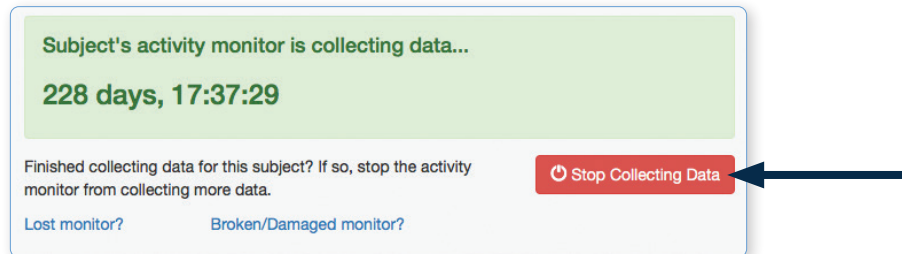

- 3 A “Data Collection Complete” tag will appear in the general tab under the “Subject Details” page.

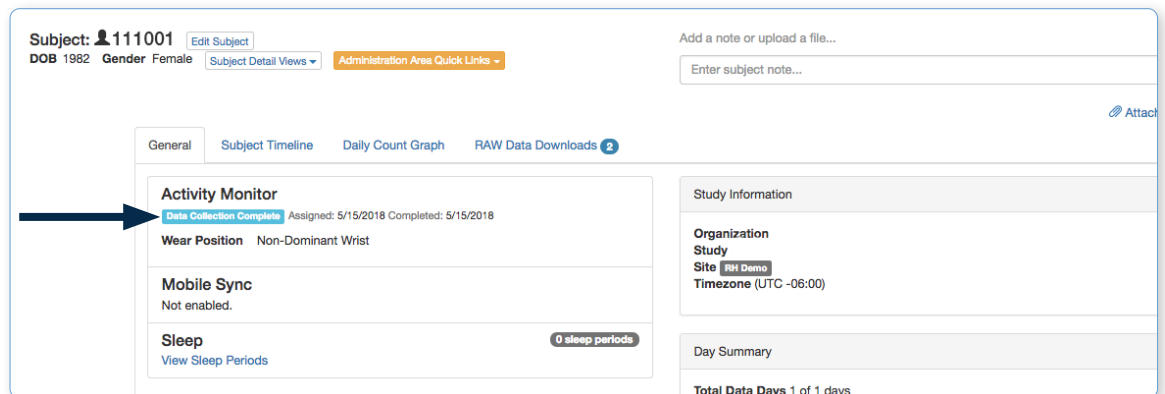


- 4 Remove the Insight Watch from the charging dock and clean with an alcohol-based solution.


Lost or Broken Activity Monitors

If an Insight Watch is lost or broken, it must be disassociated from the subject before assigning a replacement. Only one Insight Watch can be assigned to a subject at any given time.

- 1 To disassociate an Insight Watch, select the **"Subjects"** tab from the left hand menu, and select the desired subject from the list.

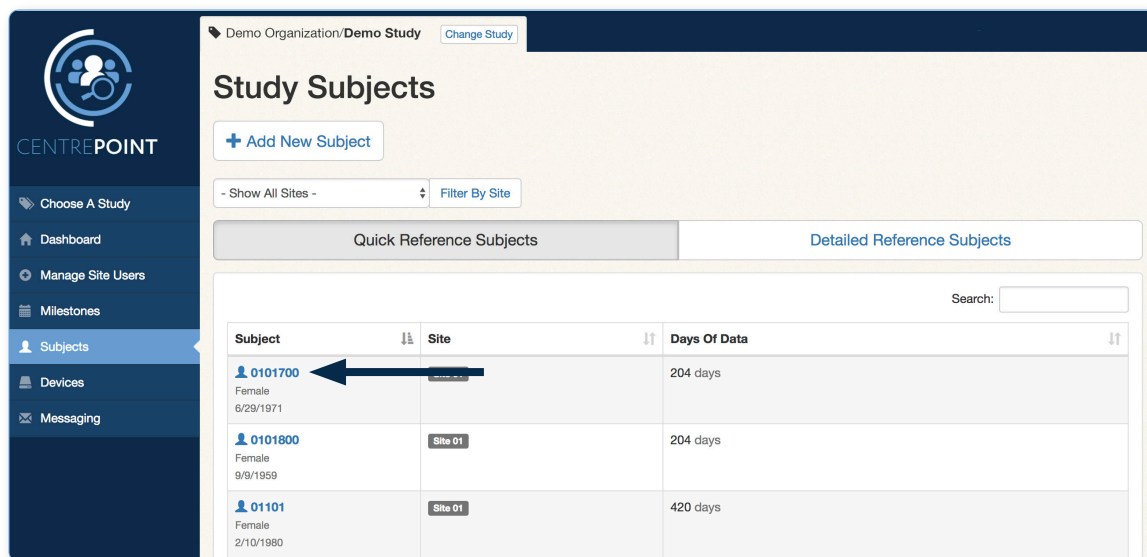

- 2 From the **"Subject Details"** page, select either **"Lost monitor"** or **"Broken/Damaged monitor."**

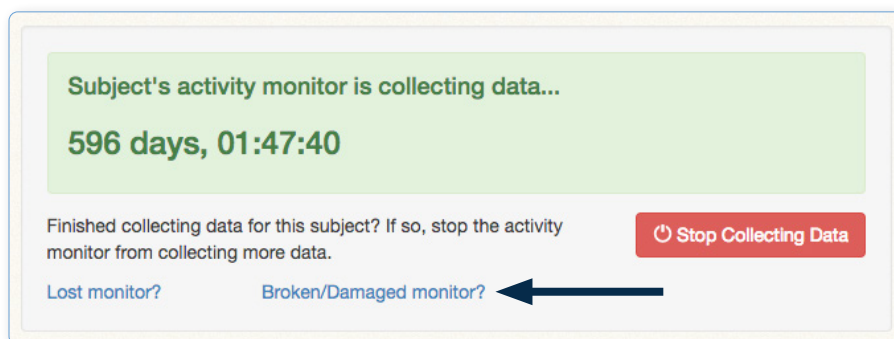

Lost or Broken Activity Monitors (Continued)

- 3 Select “Remove monitor assignment from this subject” to disassociate the Insight Watch.

Broken/Damaged Activity Monitor?

Subject	0101800	Monitor	TAS1F49160234
Study	Demo Study	Assigned	2/8/2017 1:34:02 PM
Site	Site 01	Last Activity Upload	8/30/2017 1:59:00 PM
Milestone	Visit 2	Last Raw Upload	8/30/2017 1:56:10 PM

Description of Problem

The activity monitor is broken or damaged and will not communicate with the CentrePoint software.

What Can I Do?

If the monitor assigned to this subject is known to be broken or damaged, use this page to remove the assignment of this monitor from this subject. Next, you may assign a replacement monitor to this subject by following the instructions in the [CentrePoint User Guide](#) for reassigning a monitor to an existing subject.

Is it possible to recover data from the damaged monitor?
In some cases it may be possible to recover data from a damaged monitor. Please contact [ActiGraph Support](#) or support@theactigraph.com for more information.

[Cancel and go back](#)[Remove monitor assignment from this subject](#)

- 4 To assign a replacement Insight Watch, follow the instructions on page 8 for assigning a monitor to an existing subject

Note: If you are using a CP Data Hub and disassociate an Insight Watch from a subject, the CP Data Hub will automatically disassociate as well. If a replacement Insight Watch is assigned, it must also be paired with the CP Data Hub.

Frequently Asked Questions



Is the Insight Watch waterproof?

The Insight Watch is water resistant up to 1 meter for 30 minutes according to IP57 certification. The Insight Watch can be worn during normal bathing and swimming activities.



What is the battery life of the Insight Watch?

Approximately 30 days.



When and how often does the Insight Watch need to be charged?

It is recommended that the Insight Watch be fully charged once every 3-4 weeks.



How long does it take to charge the Insight Watch?

Approximately 2-3 hours.



How should the Insight Watch be cleaned?

The Insight Watch should be cleaned with an alcohol-based solution (wipes or swabs).

Note: If an ActiGraph device is broken, malfunctioning, or does not appear to be working as intended, the subject should stop wearing it immediately and return it to the site. Site personnel should then return the device to ActiGraph.

1

Regulatory statements

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received including interference that may cause undesired operation. Changes or modifications not expressly approved by ActiGraph, LLC will void the user's authority to operate the equipment under FCC regulations.

- *FCC Part 15.107 – AC Conducted Emissions*
- *FCC Part 15.109 – Radiated Emissions*
- *FCC Part 15.207 – Modular Transmitter AC Line Conducted Emissions*
- *FCC Part 15.249 – Radiated Emission Limits of Intentional Radiators*

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- *Reorient or relocate the receiving antenna*
- *Increase the separation between the equipment and receiver*
- *Connect the equipment into an outlet on a circuit different from that to which the receiver is connected*
- *Consult the dealer or an experienced radio/TV technician for help*

2

Regulatory symbols



CAUTION: Do not simultaneously wear and charge, service, or provide any maintenance on the product.

CAUTION: Transporting or operating this product outside of the temperature range of -20°C to 40°C could lead to dangerous conditions.

CAUTION: Modification to ActiGraph products are not permitted and will void all warranties if tampered and/or modified. Do not modify the product in any way.

CAUTION: Do not swallow any part of this product. If a piece has been swallowed, contact your local poison control hotline or seek medical attention as soon as possible.

NOTICE: No precautions need to be taken in the event of changes in the performance of the product.

NOTICE: Product does not have any contraindication(s)

NOTICE: Only use ActiGraph's approved charging dock to charge the device. Do not connect the device to anything else but a computer and an ActiGraph approved charging dock.

Regulatory Information (Continued)



Emergo Europe
Prinsessegracht 20
2514 AP The Hague
The Netherlands

Authorized Representative in European Community

Indicates the authorized representative in the European Community.

AUSTRALIAN SPONSOR

Emergo Australia
Level 20, Tower II
Darling Park
201 Sussex Street
Sydney NSW 2000
Australia

ActiGraph's Australia Sponsor

Indicates ActiGraph's authorized Australian Sponsor.



Catalog Number

Indicates the manufacturer's catalogue number so that the medical device can be identified for reordering.



CE Symbol

By affixing the CE marking to a product, a manufacturer declares that the product meets all the legal requirements for CE marking and can be sold throughout the EEA. ActiGraph's products abide with the Medical Device Directive 93/42/EEC and Radio Equipment Directive 2014/53/EU.



Brazil National Telecommunications Agency (ANATEL) Homologação Number.

This equipment operates on a secondary basis, that is, not entitled to protection from harmful interference, even for stations of the same type, and may not cause interference to systems operating on a primary basis.



Consult Instructions for Use

Indicates the need for the user to consult the instructions for use.



Serial Number

Indicates the manufacturer's serial number so that a specific medical device can be identified.



Manufacturer

ActiGraph is the medical device manufacturer and is located at 49 E. Chase Street, Pensacola, Florida 32502.



Temperature Limit

Indicates the temperature limits to which the medical device can be safely exposed.



Recycle: Electronic Equipment

Indicates the medical device should not be disposed of in the trash. Contact ActiGraph Customer Service regarding the disposal of these products.



Lithium Ion (EU)

Indicates the lithium ion battery within the device should not be disposed of in the trash. Contact ActiGraph Customer Service regarding the disposal of these batteries.



ActiGraphCorp.com

Please contact our Customer Support Team with any questions or for additional information about operating the ActiGraph CentrePoint Insight Watch.

49 East Chase Street. Pensacola, FL 32502

support@actigraphcorp.com | 850.332.7900 | +1 (877) 497-6996 (toll free)

©2019 ActiGraph, LLC. All rights reserved.

U.S. federal law restricts this device to sale by or on the order of a physician.